


H4™ Pitts BROAD ARCHWIRE SEQUENCE

| WEEKS | NON-EXTRACTION | | FLIPPED NON-EXTRACTION | | EXTRACTION | | | |
|---|--|--------------|--|--------------|---|--------------|------------------------------|--|
| 5 | .014 Thermal Activated (6 weeks) | | | | | | | |
| 10 | | | | | | | | |
| 15 | | | | | | | | |
| 20 | .018 x .018 Ultra Soft Thermal Activated (6 weeks) | | | | | | | |
| 25 | | | | | | | | |
| 30 | | | | | | | | |
| PAN-REPO ADJUST CASE MANAGEMENT (PRACM) | | | | | | | | |
| | Rotations | No Rotations | Rotations | No Rotations | Rotations | No Rotations | | |
| 35 | RESET BRACKETS OR PICK UP SECOND MOLARS | | .020 x .020 Thermal Activated | | | | | |
| 40 | .019 x .025 Beta Titanium | | | | FOR ENMASSE SPACE CLOSURE .018 x .025 Stainless Steel | | | |
| 45 | | | | | | | .020 x .020 Beta Titanium | |
| 50+ | | | | | .020 x .020 Beta Titanium | | .020 x .020 Beta Titanium | |
| TO DEACTIVATE THE APPLIANCE WHEN IDEAL AXIAL INCLINATION ATTAINED OR BIOLOGICAL LIMIT IS REACHED | | | | | | | | |

Courtesy of: Dr. Tom Pitts and Dr. Duncan Brown

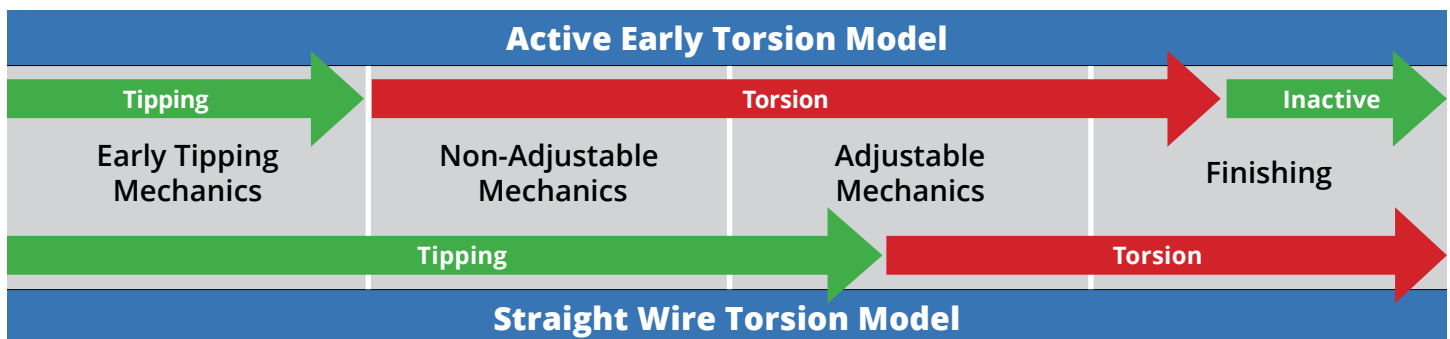

Courtesy of: Dr. Tom Pitts and Dr. Duncan Brown